

Server Update

For eRAD RIS
Version 1.0
Build 41
Update 41.003

Table of Contents

1. New Features	3
Reporting - Limit or convert special characters in cue names, macros and templates	3
Reporting - Radiologist edit text color.....	4
Reporting - field tooltips staying on top of other screens after displaying	5
2. Defects Resolved	6
3. Code Stream.....	7
4. Package Contents	8
5. Intended Audience	9
6. Who is affected	9
7. eRAD RIS Release Version Numbers	10
8. Installing	11
9. Known Limitations	11

1. New Features

Reporting - Limit or convert special characters in cue names, macros and templates

Conversion of special characters in cue names, macros and templates.

In a template if you have a field with a value in the cue like this "Technique :"

When you save the template and reopen it, the : is converted as follows "Technique %colon%".

This feature does the following, while dictating an exam and you say the words Technique colon, your cursor will be brought to the field that has a value of "Technique %colon%" in the cue.

NOTE: The user preference "Aggressive field navigation" must be checked for this feature to work, after checking this preference, a restart of RIS is required.

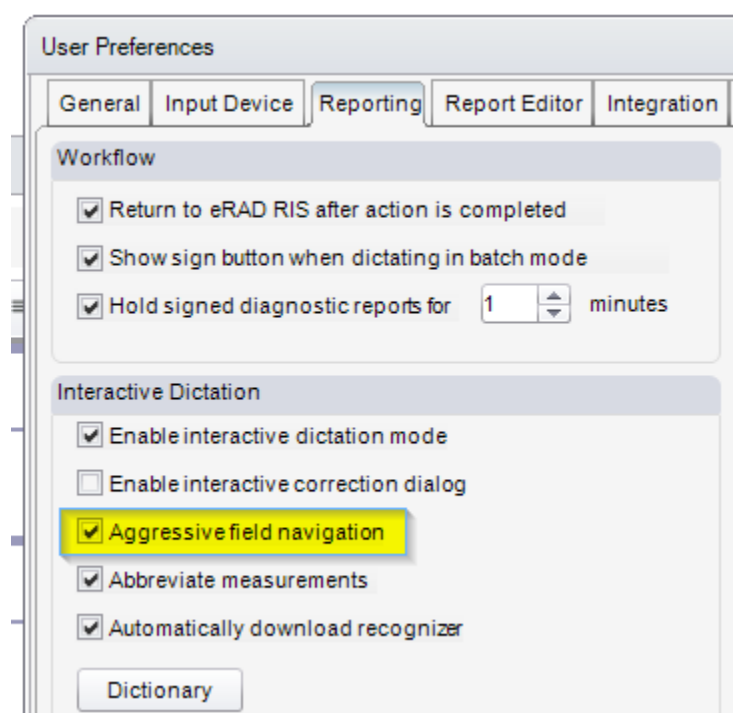


Figure 1.1 – User preference “Aggressive field navigation”

NOTE: the user preference “Aggressive field navigation” will need to be selected for “No Wait Navigate”. We had issues in the past when dictating in Field 1 and before the text was inserted the command for next field was issued. The text that had not yet been inserted in Field 1 would appear in Field 2. This has been addressed, but as mentioned above the Aggressive field navigation user preference needs to be enabled.

Reporting - Radiologist edit text color

Within the reporting window, the user can open the context menu (right click) on the editor. This will present the user with a list of options. Selecting the “Editor Appearance” option will open the “Editor Appearance” window. A new option has been added.

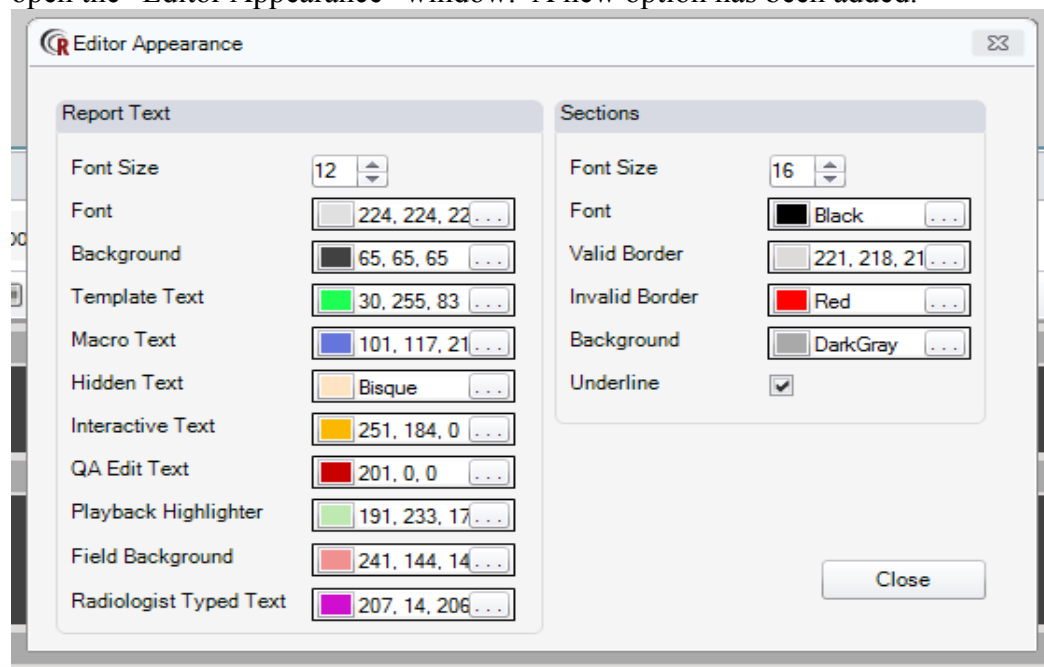


Figure 1.2 – Editor Appearance window with new Radiologist Typed Text option

This new option controls the color of the manually entered/edited text by the radiologist

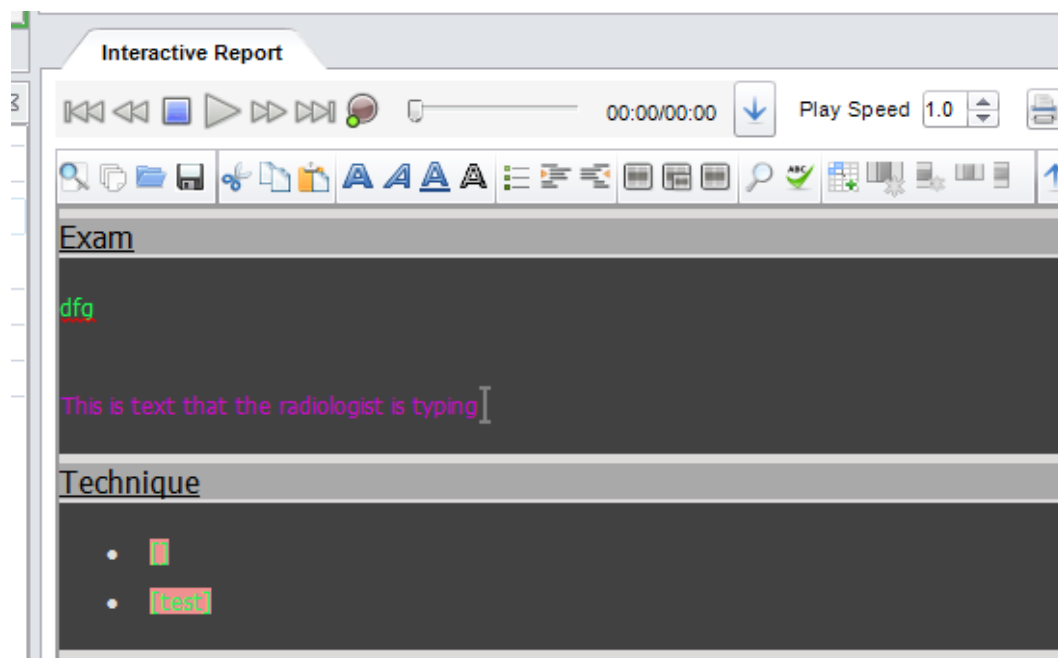


Figure 1.3 – Editor displaying the text color when manually entered by the radiologist.

Reporting - field tooltips staying on top of other screens after displaying

Templates used in the reporting window can have navigational fields pre-defined in the template. Each of these fields has a name associated with it. If the mouse pointer hovers over a field it will display the field name in a tooltip.

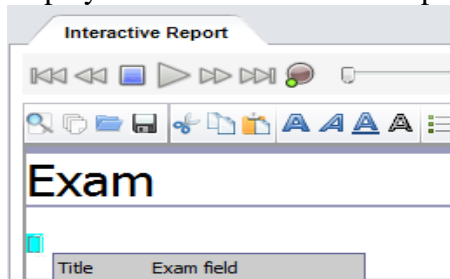


Figure 1.4 – Reporting field tooltip.

A new user preference on the Report Editor tab called “Show report field tooltips” will enable/disable the tooltip from displaying.

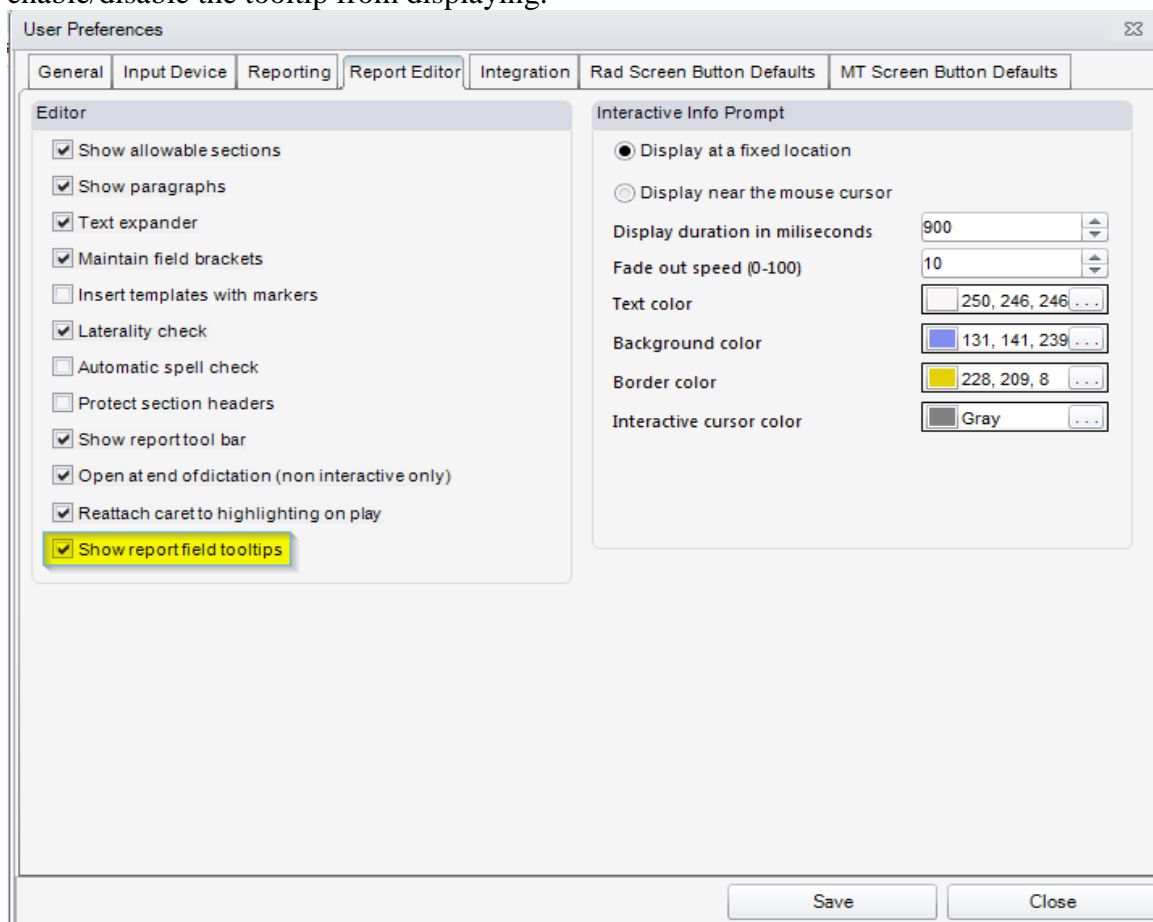


Figure 1.5 – New User Preference “Show report field tooltips”

RADNET, Inc. – eRAD RIS Server Update

2. Defects Resolved


Redmine #	Subject	Priority
5251	Document Distribution - EMR Sender functionality broke in B41	Immediate
5163	Losing default value for fields when saving template/macro	Immediate
5212	Scanner loads previous patients scans in preview after scanning (RRR)	Urgent
5184	Arrive individual exams error (RRR)	Urgent
5035	Wrong document model is loading in template editor	Urgent
5140	Editing fields in template editor sometimes throws an exception	High
5036	Reporting - ability to navigate cued fields without pausing	High
5034	Reporting - 7.82+ Numbering paragraphs creates an unwanted list	High
5032	Reporting - playback is not in order of dictated	High
5030	Sections from incorrect practice loading into template editor	High
5027	Using next field button while dictating (Fisher site visit)	High
5167	Aggressive field cue navigation fails after editing a cue on a template.	Normal
5149	Reporting - speech mics 3500 and 3600 showing incorrectly in Help>About	Normal
5049	Reporting - New recognizer for 7.88	Normal
5047	Next field voice command, showing field popup for one field too far.	Normal
5046	Reporting - ability to configure the completeTimeout setting for interactive dictation	Normal
5045	Reporting - ability to configure the finalTimeout setting for interactive dictation	Normal
5044	Reporting - ability to configure the speedVsAccuracy setting for interactive dictation	Normal
5043	Reporting - improper capitalization of measurements outside of a field	Normal
5042	Reporting - structured lists splitting when dictating next item in the middle of an existing item	Normal
5041	Reporting - ability to dictate 'next item' without pausing	Normal
5040	Reporting - integrate M*MOdal 7.88 release into 41.1 code stream.	Normal
5039	Reporting - negative two views rendering	Normal
5038	Reporting - Limit or convert special characters in cue names, macros and templates	Normal
5037	Reporting - ability to dictate negative or minus numbers	Normal
5033	Reporting - interactive dictation of the word "through" appearing as "-" not ideal	Normal
5031	Reporting - AudioReplacementEnabled on causes "undo that" to malfunction	Normal

RADNET, Inc. – eRAD RIS Server Update

5028	Improve performance of "next field" (CRA site visit)	Normal
5023	Reporting - Not a field exception causing application hang	Normal
5022	MModal - select command will not select first or last word in a field	Normal
5020	Reporting - interactive number lists are not consistent	Normal
5019	Reporting - field tooltips staying on top of other screens after displaying	Normal
4801	NewCrop UPIN problem	Normal
5029	Unable to dictate "+" into a report (CRA site visit)	Low

3. Code Stream

eRad RIS– code merge

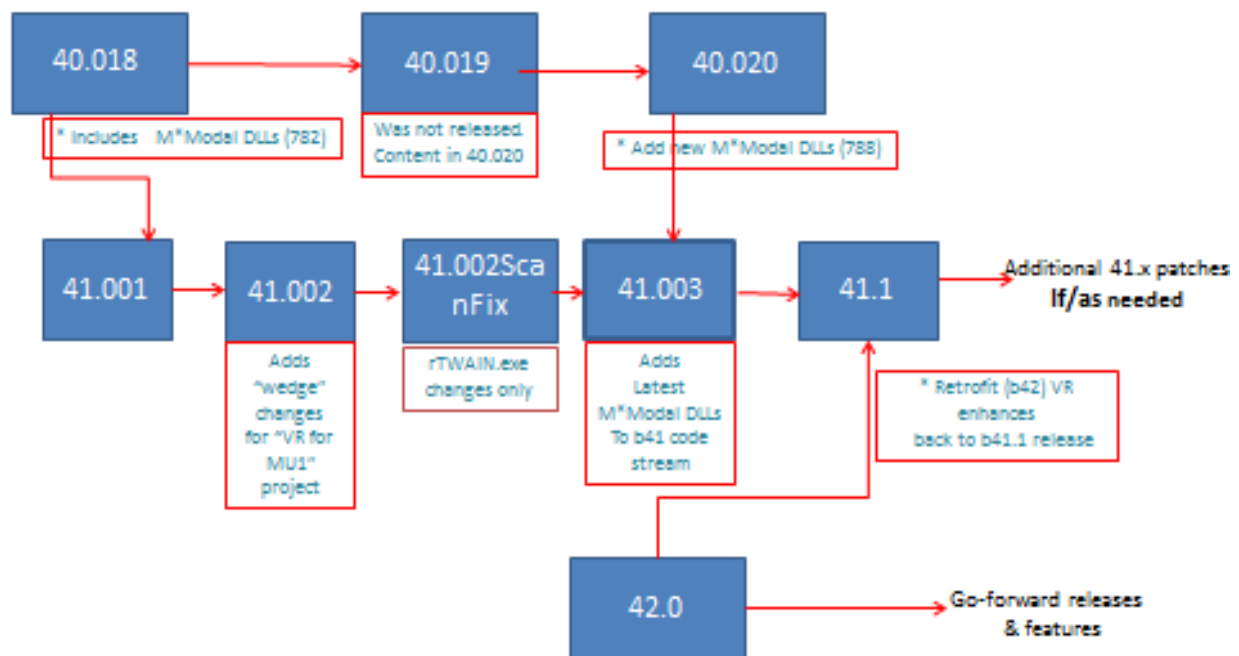

4. Package Contents

Files included for this server update:

- 1) GUI ZIP Updates folder contains
 - a) rRIS_1.3.41.20171 (3GB).zip

- 2) RIS Service Updates folder contains the following files.
 - a) CDSJobManager.dll
 - b) CDSJobManager.pdb
 - c) CommonLib.dll
 - d) CommonLib.pdb
 - e) DDEngine.dll
 - f) DDEngine.pdb
 - g) rRISServices.exe
 - h) rRISServices.pdb
 - i) rRISServicesx64.exe
 - j) ScheduledTaskManager.dll
 - k) ScheduledTaskManager.pdb
 - l) SignatureApprovalService.exe
 - m) SignatureApprovalService.pdb
 - n) WCFService.dll
 - o) WCFService.pdb
 - p) WebApi.dll
 - q) WebApi.pdb
 - r) Wedge.exe
 - s) Wedge.pdb

- 3) Recognizer folder contains:
 - a) 7~Recognizer.MtiRecognitionUS.Radiology_Interactive_v13_v5_401.12965.xc
 - b) 7~Recognizer.MtiRecognitionUS.Radiology_Interactive_v13_v5_401.12965-descriptor.xc

5. Intended Audience

This document is created by the RIS Development team for the RadNet RIS Implementation team.

6. Who is affected

Build 41 installs. This server update must be applied directly to the build 41 core install. Server updates are not cumulative. Any previous server updates must be applied before this one.

RADNET, Inc. – eRAD RIS Server Update

7. eRAD RIS Release Version Numbers

Build	Patch	UI Version	Core Version	WS Version	DB Version	Notes
40	-	1.3.40.16628	1.3.40.16628	1.3.40.16628	1.3.40.16628	Full version release
40	1	1.3.40.16849	1.3.40.16659	1.3.40.16744	1.3.40.16727	GUI.zip, Web service, database updates and SQL Management reporting
40	2	1.3.40.16868	1.3.40.16659	1.3.40.16868	1.3.40.16727	GUI.zip and Web Service
40	3	1.3.40.16868	1.3.40.16659	1.3.40.16868	1.3.40.16727	DB updated only
40	4	1.3.40.16967	1.3.40.16659	1.3.40.16932	1.3.40.16961	GUI.zip, Web service, database updates
40	5	1.3.40.17087	1.3.40.16659	1.3.40.17069	1.3.40.17087	GUI.zip, Web service, database updates
40	6	1.3.40.17146	1.3.40.16659	1.3.40.17146	1.3.40.17146	GUI.zip, Web service, database updates
40	7	1.3.40.17268	1.3.40.16659	1.3.40.17268	1.3.40.17259	GUI.zip, Web service, database updates
40	8	1.3.40.17364	1.3.40.16659	1.3.40.17362	1.3.40.17343	GUI.zip, Web service, database updates
40	9	1.3.40.17370	1.3.40.16659	1.3.40.17370	1.3.40.17370	GUI.zip, Web service, database updates
40	10	1.3.40.17435	1.3.40.16659	1.3.40.17435	1.3.40.17370	64 BIT web service. GUI.zip, Webservice.
40	11	1.3.40.17435	1.3.40.16659	1.3.40.17525	1.3.40.17370	Web service update
40	12	1.3.40.17685	1.3.40.16659	1.3.40.17685	1.3.40.17648	GUI.zip, Web service, database updates
40	13	1.3.40.17776	1.3.40.16659	1.3.40.17776	1.3.40.17648	GUI.zip, Web service
40	14	1.3.40.18063(3GB)	1.3.40.16659	1.3.40.18063	1.3.40.18063	GUI.zip, Web service, database updates
40	14.1	1.3.40.18063(3GB)	1.3.40.16659	1.3.40.18107	1.3.40.18063	Web service
40	15	1.3.40.18127(3GB)	1.3.40.16659	1.3.40.18127	1.3.40.18127	GUI.zip, Web service, database updates
40	16	1.1.3.40.18585 (3GB)	1.3.40.16659	1.3.40.18585	1.3.40.18646	GUI.zip, Web service, database updates
40	16.1	1.3.40.18705 (3GB)	1.3.40.16659	1.3.40.18585	1.3.40.18646	GUI.zip
40	17	1.3.40.18780 (3GB)	1.3.40.18756	1.3.40.18756	1.3.40.18646	Core install required for Radiologist User Group. GUI.zip, Core setup, and Web Service updates
40	18	1.3.40.19253(3GB)	1.3.40.18756	1.3.40.18756	1.3.40.19253	GUI.zip, Web Service and database updates
40	19	1.3.40.19909(3GB)	1.3.40.18756	1.3.40.19878	1.3.40.19253	GUI.zip, Web Service
40	20	1.3.40.20086(3GB)	1.3.40.18756	1.3.40.19878	1.3.40.19253	Gui.zip and Recognizer
41	-	1.3.41.19135(3GB)	1.3.41.19135	1.3.41.19135	1.3.41.19135	Full version release
41	1	1.3.41.19356(3GB)	1.3.41.19135	1.3.41.19346	1.3.41.19346	GUI.zip, Web Service and database updates
41	2	1.3.41.19620(3GB)	1.3.41.19135	1.3.41.19620	1.3.41.19346	GUI.zip and Web Service
41	3	1.3.41.20171(3GB)	1.3.41.19135	1.3.41.20171	1.3.41.19346	GUI.zip , Web Service and Recognizer

8. Installing

- 1) Copy and replace the current rRIS...zip file with the rRIS_1.3.41.20171.zip file provided with this release.

Be sure to make a backup of the ris.exe.config file contained within the current .zip folder.

- 2) Stop each service listed below. Copy the files from the RIS Service Updates folder to the directory for each rRISService. Be sure to make backup of all files, but especially the rRISServices.exe.config and rRISServices.vshost.exe.config files.

rRISService directories include:

- a) *rRISService*
- b) *CDS rRISService*
- c) *Wedge rRISService*
- d) *DDE rRISService*

Restart each service.

- 3) Copy the Recognizer files to your local machine to the following folder.

Note: this could be a hidden folder.

C:\ProgramData\RadNet\Recognizer\2.16.840.1.113883.3.21\

9. Known Limitations

Known limitations listed in this document may not have been introduced in 41.003, but have been discovered during testing of 41.003.

Redmine #	Subject	Info
5245	Reporting - Spine navigational cues are not recognized	<p>My cues were configured as</p> <ol style="list-style-type: none">1. C two dash three colon2. C three dash four colon3. C four dash five colon <p>I changed them to be</p> <ol style="list-style-type: none">1. C two colon2. C three colon3. C four colon <p>The change in cues names worked much better</p>
5248	Reporting - Edit template from reporting window error	<p>This was only occurring with one of my templates. If this occurs try "refresh configuration". If that does not work</p>

RADNET, Inc. – eRAD RIS Server Update

		open the same template in the TemplateMacro editor and save.
5249	Reporting - Shift Tab in the template editor or reporting editor will return error	
5219	Reporting - Function key action DeleteWordLeftOfCursor will delete fields and section titles	If you continue to select the function key after the text has been removed it will jump to the section title and can also delete it as well
5220	Reporting - Interactive Info Prompt can show fields that have been jumped over when using Next Section command	If in a section that has multiple field with default text, issue the command Next Section, the interactive info prompt will display the text from the fields that were skipped.
5254	Reporting - using the undo feature, the cursor is not following the text that is being removed.	Seems to occur when the “undo” feature is used multiple times concurrently.
5260	Skip and Continue only loads images every other study	Two workarounds exist for this. 1. Turn off (unselect) user preference “Reuse reporting screen” 2. If using “Reuse reporting screen” select the Load images icon in the reporting window.