



eRAD RIS

RELEASE ANNOUNCEMENT

Rapid Release 4.2026.002.2.1

UPDATED MARCH 17, 2026

TABLE OF CONTENTS

Summary 2

 Release Announcement2

 Resolved Issues Error! Bookmark not defined.

Version Details 2

 Code Stream4

PUBLICATION HISTORY

Revision	Author	Description
March 17, 2026	Kevin Brooks	<ul style="list-style-type: none"> Commercial release.
March 18, 2026	Kevin Brooks	<ul style="list-style-type: none"> Corrected release version from 4.2026.002.3.1 to ase 4.2026.002.2.1

SUMMARY

Release Announcement

This Rapid Release release of ERAD RIS 4.2026.002.2.1 introduces an enhancement to Patient Portal to strip unexpected delimiters from migrated Provider Portal user account data.

New Features

This release introduces the following features and enhancements:

Category	Redmine #	Subject	Description
Patient Portal	38218	Account migration from Akumin Millinger RF portal to erad RF portal decryption data parsing issue	This enhancement to Patient Portal introduces a mechanism to strip unexpected delimiters from migrated Provider Portal user account data.

Refer to the **FEATURE DETAILS** section below for configuration and usage information.

FEATURE DETAILS

Patient Portal

Feature #38218 - Portal User Account Migration Data Parsing Issue

Summary

This enhancement to Patient Portal introduces a mechanism to strip unexpected delimiters from migrated Provider Portal user account data.

Background

While migrating Provider Portal account data Millinger RF Portal to eRAD Provider Portal, unexpected characters in the data (specifically the ^{AMP} symbol) resulted in decryption data parsing errors.

Feature Description

With this change, two new configuration settings have been introduced, identifying a target and replacement text for use during data migration.

Service Team Upgrade Instructions

The Service Team must complete the following actions to deploy this feature:

Portals

Changes to applicationsettings.config

- Enable this feature by configuring the two configuration values.

Setting	Default Value	Purpose
AccountMigrationRISIDPlaceholder	^AMP^	Migration data source (problem) text.
AccountMigrationRISIDEscapedChar	&	Migration data replacement (fixed) text.


Configuration Instructions

Service Team assistance is required to enable this feature.

VERSION DETAILS

Code Stream

The following source code branches have been merged into this release:

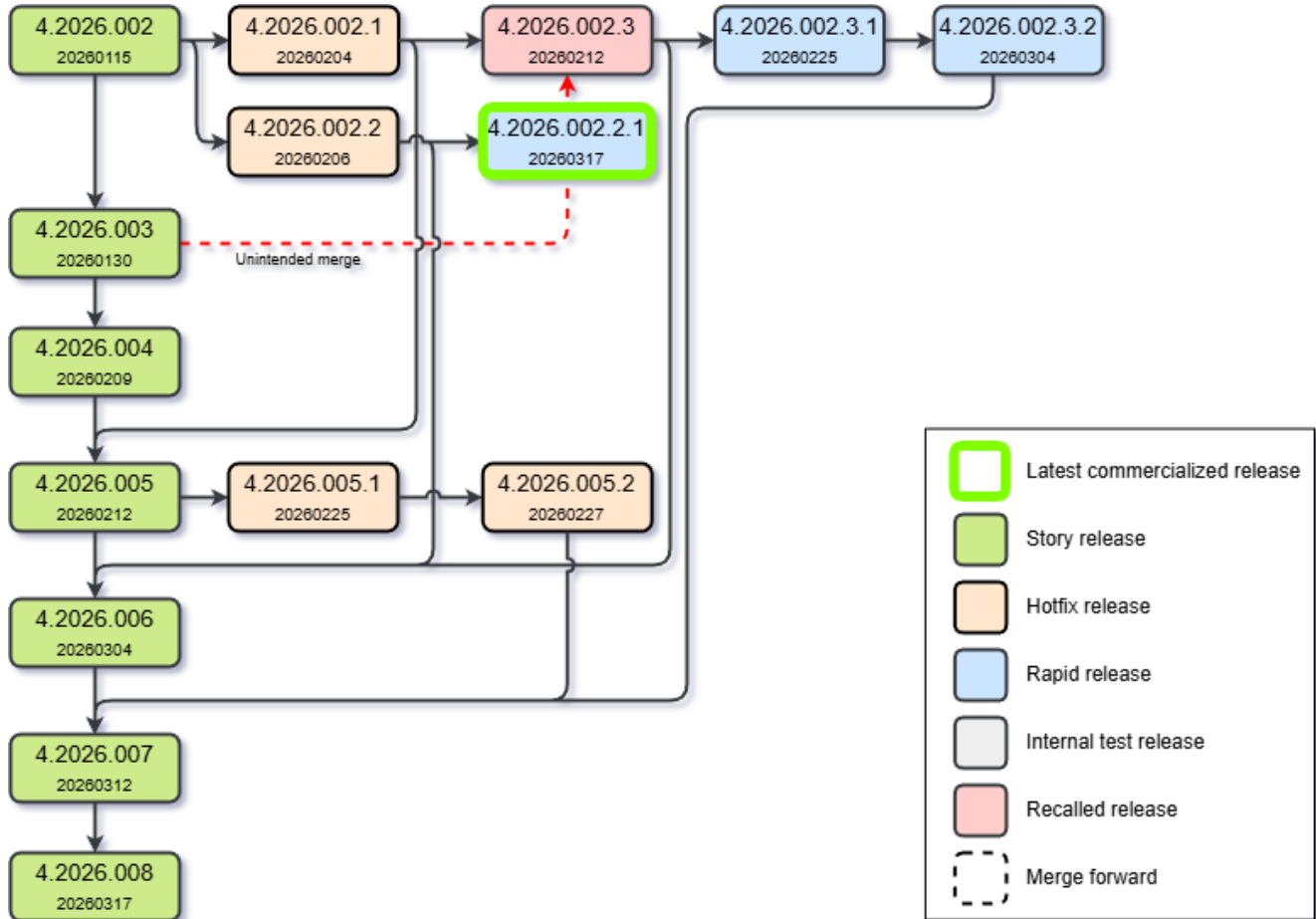

WARNING: The functionality of this release has not yet been merged forward into any other release.

Please carefully review the code stream diagram of any future releases before deploying to ensure that the functionality of this release will not be lost when upgrading.