


For some types of exams, the radiologist may want to enter a finding (positive or negative) for which a set of follow-up recommendations is available. Similar to the Mammo Tracking workflow, this allows the staff to work cases requiring follow-up from a designated worklist under the Radiologist menu, called the Finding Follow Up worklist.

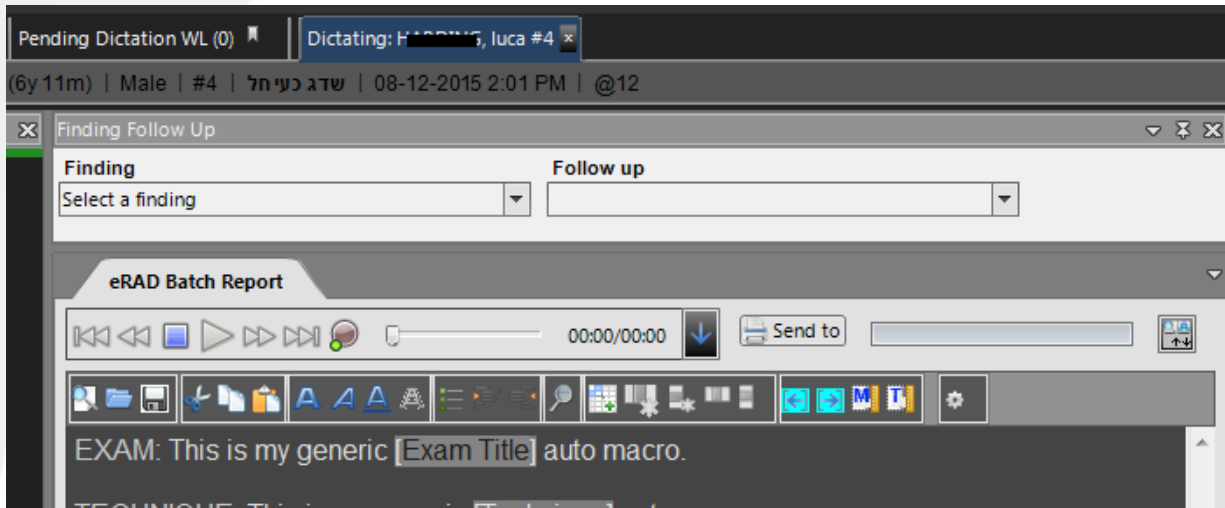

Note that if a finding is required for that exam, the radiologist will not be able to sign the report until the Finding drop-down is populated. However, the radiologist can send the report to transcription, and the editor can populate the value based on the dictation.

The options available in the finding drop-down are populated from entries in the Finding lookup table (under Configuration).

| Finding Code | Positive Finding Flag | Display Order | Description |
|-----------------------------|-----------------------|---------------|--|
| Click here to add a new row | | | |
| Rad | N | 101 | Radiologist |
| Tech | N | 100 | Technologist |
| Brain lesion | Y | 98 | Brain lesion(s) |
| Comparison | N | 99 | Comparison with Priors Required |
| LUNG 4X | Y | 25 | Suspicious |
| LUNG 4Xc | Y | 27 | Suspicious, prior hx of lung cancer |
| LUNG 4Xs | Y | 26 | Suspicious with other significant finding |
| LUNG 4Xsc | Y | 28 | Suspicious with other significant finding, prior hx of lung cancer |
| LUNG 4Ac | Y | 19 | Suspicious, prior hx of lung cancer |
| LUNG 4As | Y | 18 | Suspicious with other significant finding |

The follow-up choices can have a configurable time frame associated with each option.

Finding Follow Up (1 of 3) x

| Site | Follow Up Date | Minutes Until Follow Up | Finding | Follow Up |
|-------|--------------------|-------------------------|--|----------------------|
| | | | | StartsWith: BX |
| Cedar | 03-16-2016 9:47 PM | -57 | LUNG 4Xs - Suspicious with other significant finding | BX - Tissue sampling |

Any exams that require follow up will be placed onto the worklist for a user to review and resolve.

Finding Follow Up (6) x

| Flags | Status | Procedures | Referring | Finding | Follow Up | Time Until Follow Up | Minutes Until Follow Up | Follo |
|--------|--------|---|--------------------|--|--|----------------------|-------------------------|-------|
| Signed | | CT Smoker's Lung Screening [71250S] - Chest | Tham, Dale | LUNG 3s - Probably benign with other significant finding | 6-M - Six month LDCT | -162d 19h 50m | -234471 | 12-2 |
| Signed | | CT Screen Chest Wo [CT39] - Chest | Washington, George | LUNG 1c - Negative, prior hx of lung cancer | ANNUAL - Annual screening with LDCT in 12 months | -82d 7h 17m | -118518 | 03-1 |
| Signed | | CT Smoker's Lung Screening [71250S] - Chest | Test, Henry | LUNG 4Xs - Suspicious with other significant finding | BX - Tissue sampling | -448d 7h 56m | -645596 | 03-1 |
| Signed | | CT Smoker's Lung Screening [71250S] - Chest | Waite, Andrew R | LUNG 3 - Probably benign | 6-M - Six month LDCT | -288d 7h 28m | -386369 | 09-1 |
| Signed | | CT Smoker's Lung Screening [71250S] - Chest | Warren, Spencer | LUNG 1 - Negative | ANNUAL - Annual screening with LDCT in 12 months | 23d 5h 35m | 33455 | 07-0 |
| Signed | | CT Smoker's Lung Screening [71250S] - Chest | Franklin, Ben | LUNG 0 - Incomplete | COMP - Comparison with prior chest CT | -342d 18h 19m | -493579 | 06-3 |

From the Finding Follow Up worklist, manage a follow-up by right-clicking the row and selecting Remove Follow Up Required Flag. Or open the Follow Up screen shown below and select the Resolved checkbox at the bottom.

eRAD RIS (Development as Denis test Doiron [denis])

File Administration Messaging Scheduling Front Desk Technologist Radiologist Mammo Editor PACS Help

All Pending Dictation WL (2) Pending Dictation WL (0) Finding Follow Up (1) Follow Up: HARRIS, luca #4

Patient Order Contact Log

Contact Log

| Notes | Date | Contact by | Type |
|-------|------|------------|------|
| | | | |

Add View

Patient Notes

08-14-15 denis -

Previous Patient Notes

Resolved