


## Pathology Results

All signed cases with BI-RADS of 4 or 5 (biopsy recommendations) will be listed on the ACR Category 4/5 Follow up WL under the Mammo Menu.

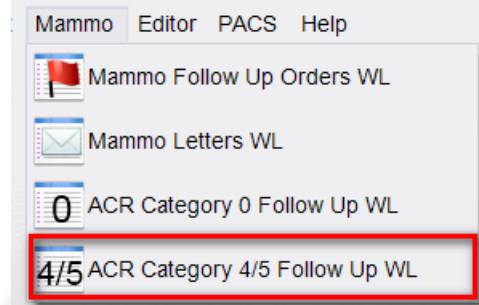

BIRAD 4/5 patients will be listed here

ACR Cat. 4/5 Follow-Up WL (2) x

(use location filter) BI-RADS 4/5 DEMO BI-RADS 4/5


Flags	Status	Procedures	Scheduled Date	Referring	First Name
	Signed	MA Digital Mammo Screen Bilat [G0202] - Bilat...	06-03-2007 5:09 PM	Wilson, John	Allison
	Signed	MA Digital Mammo Diag Bilat [G0204] - Breast -...	03-11-2016 10:38 AM	Washington, George	Hilary

Click the Follow-Up button or double click the patient to:

- Track attempts to obtain pathology results.
- Add pathology data once received.
- Please note that recommendations made in a Legacy System must be tracked in the system from which the recommendation originated. The pathology report can be scanned into eRAD RIS for reference and if the radiologist wishes to dictate an addendum to the report, this can be done in eRAD as well.

Follow Up

The Contact Log should be used to record attempts to obtain pathology results for biopsies performed at an outside facility. There is a Contact Type called Pathology Requested.

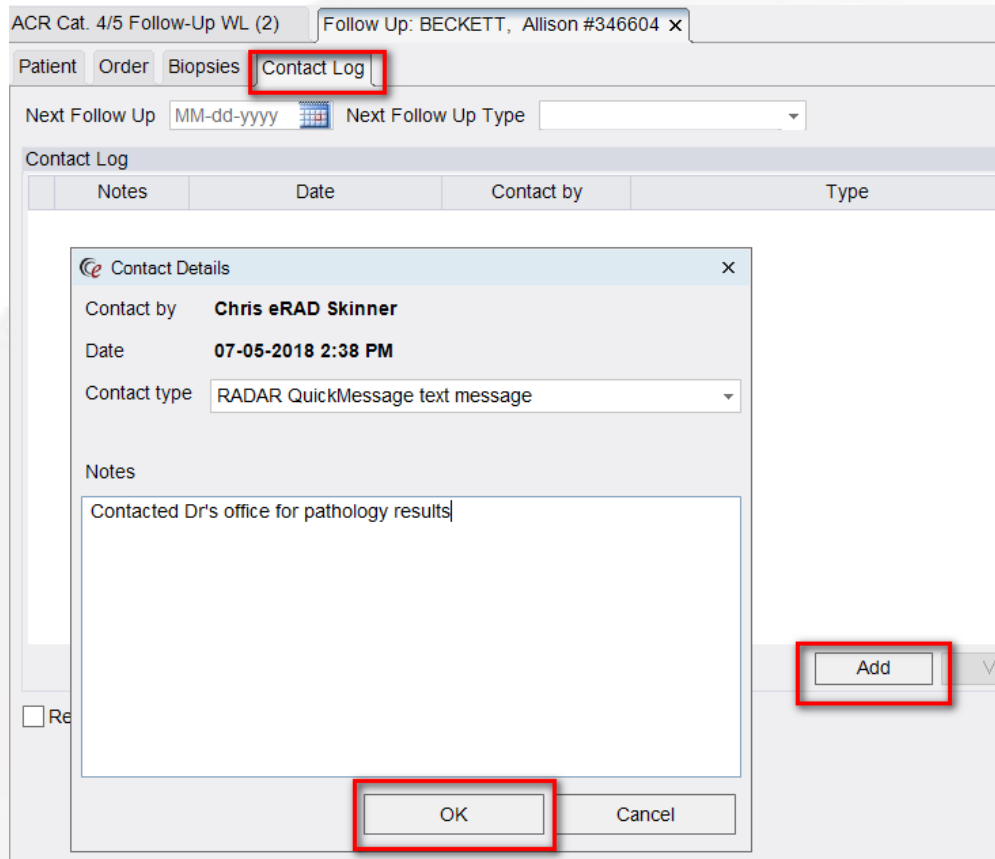

ACR Cat: 4/5 Follow-Up WL (2)    Follow Up: BECKETT, Allison #346604 x

Patient    Order    Biopsies    **Contact Log**

Next Follow Up: MM-dd-yyyy    Next Follow Up Type: [v]

Notes	Date	Contact by	Type
<div style="border: 1px solid gray; padding: 5px;"> <p><b>Contact Details</b> [x]</p> <p>Contact by: <b>Chris eRAD Skinner</b></p> <p>Date: <b>07-05-2018 2:38 PM</b></p> <p>Contact type: RADAR QuickMessage text message [v]</p> <p>Notes: Contacted Dr's office for pathology results</p> </div>			


Resolved

**Add**    View

**OK**    Cancel

If the case is to be resolved without a biopsy being performed, due to the patient or referring physician choosing an alternate course of action or if you have not been able to obtain pathology results after your designated number of attempts, enter an appropriate Contact Log entry and then check the Resolved checkbox on the Contact Log tab. Choose the corresponding reason on the pop-up window.

Resolved

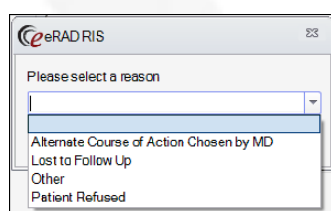

Please select a reason

- Alternate Course of Action Chosen by MD
- Lost to Follow Up
- Other
- Patient Refused

Follow Up: BECKETT, Allison #346604 x


Patient Order **Biopsies** Contact Log

Mammo Biopsies

Date performed	Technique	Laterality


Add Edit Remove

To add a biopsy to the record, click the Add button and enter the following information:

- ☞ Performing Surgeon -- To select the radiologist or other physician who performed the biopsy, we must have that doctor marked as a surgeon in eRAD. All of the surgeons are listed in the drop-down menu alphabetically.
- ☞ Biopsy Location -- Enter the name of the facility where the biopsy was performed.
- ☞ Biopsy Technique -- Click the drop-down to select a technique from the list.
- ☞ Date Performed -- Enter the date of the procedure or click the calendar button  to select a date.
- ☞ Lesion Classification -- Choose Benign, High Risk Benign, or Malignant to classify the lesion.
- ☞ Breast Side -- Choose which breast was biopsied.
- ☞ Performed -- Indicate where the biopsy was performed.
- ☞ Additional Surgery -- Indicate whether further pathology from additional surgery or a follow-up biopsy is required *for a diagnosis*. This is usually only selected if the pathology results are not concordant with the imaging findings or the sample is insufficient and repeat biopsy is recommended. However, some sites keep the patient on the list to obtain Surgical Pathology results in addition to the MQSA-required pathology. The user must choose one of these options:
  - Required - Patient will remain on the BI-RADS 4/5 Worklist until additional pathology is entered
  - Not Required - Patient will drop from the BI-RADS 4/5 Worklist and the Follow Up option below becomes active

For Follow Up, choose either Required or Not Required.

- ☞ If Required is selected, you will see the following pop-up window:

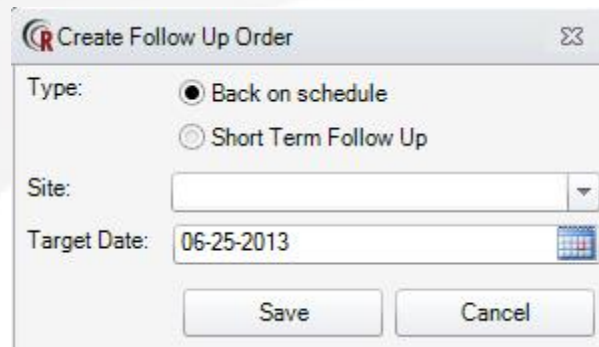


- ☞ Depending on the results of the exam and the radiologist's recommendation, you can choose either Back on Schedule for the patient to receive a Reminder Letter for their annual exam *OR* Short Term Follow Up to generate a Reminder Letter for a diagnostic mammogram or breast ultrasound to be performed at a shorter interval. Either way, be sure that the correct Target Date is chosen in order to generate the Reminder Letter appropriately.
- ☞ If Not required is selected for an ACR Category 4/5 exam, the follow-up auto-resolves from the Worklist and the patient will not receive a reminder letter for a future breast imaging exam.
  - There is a Notes field that can be used for additional information.

Technical: Fill in any applicable information for the fields on the right side of the screen if desired. The only *operationally* required field in the Technical section is "Pathology Codes;" see below for details.


- ☞ Nipple Involved -- Indicate whether the biopsy involved a nipple.
- ☞ Tumor Size -- Enter the size of tumor, in millimeters, in this field.
- ☞ Nodes Removed -- Enter the number of nodes removed.
- ☞ Nodes Positive -- Enter the number of nodes that were positive.
- ☞ Estrogen Receptor -- Indicate whether the estrogen receptor was positive or negative.
- ☞ Progesterone Receptor -- Indicate whether the progesterone receptor was positive or negative.
- ☞ Stage -- Enter the clinical stage of any malignancy.
- ☞ S Phase -- Enter the S phase of the cells.

- ☞ Margin Status -- Indicate whether the margin status is involved or uninvolved.
- ☞ **Pathology Codes** -- Enter any pathology codes that correspond to the biopsy. Start typing the pathology description to find a match, or click in the field and hit space bar to see a list of all possible choices. More than one selection can be made.
- ☞ Differentiation -- Enter the differentiation: well, moderate or poor.

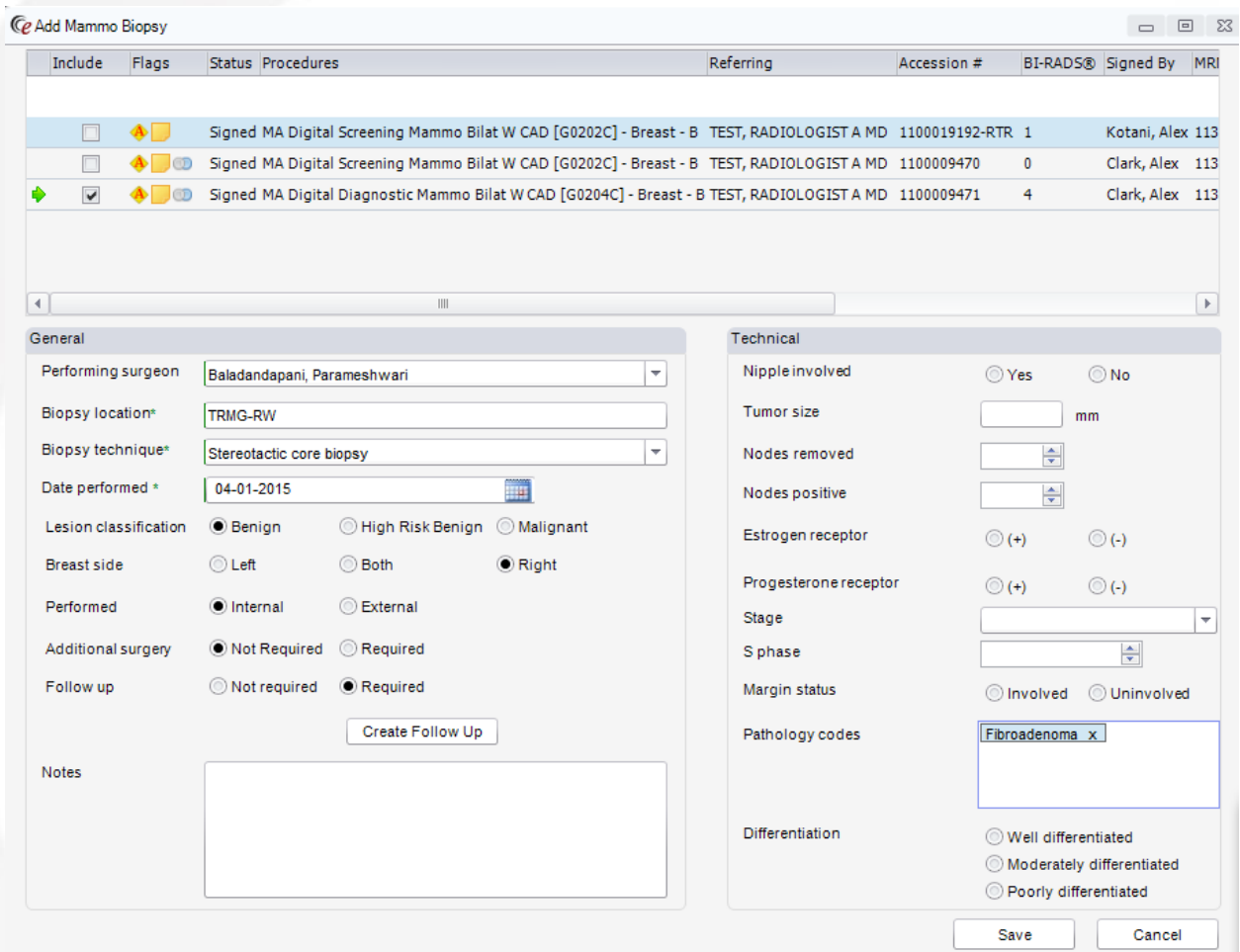

- ☞ Before saving, please note that you have the option to copy the biopsy information to other BI-RADS enabled exams, such as mammography, breast ultrasound, or breast MRI. If you do so, the information will be added to the other exam's biopsy grid. However, if the other exam is on a worklist (BI-RADS 0 or BI-RADS 4/5) it will not automatically resolve the entry, in case additional steps are required.

- In the example illustrated above, you might consider checking the box for the BI-RADS 0 screening that occurred before the BI-RADS 4 recommendation, since the BI-RADS 0 ultimately lead to the biopsy recommendation.
- ☞ Click the Save button to add the biopsy information. When saved, the biopsy is added as a row to the Mammo Biopsies grid. Click the Edit button to modify it, or click the Remove button to delete it. You may add a second biopsy by clicking Add.

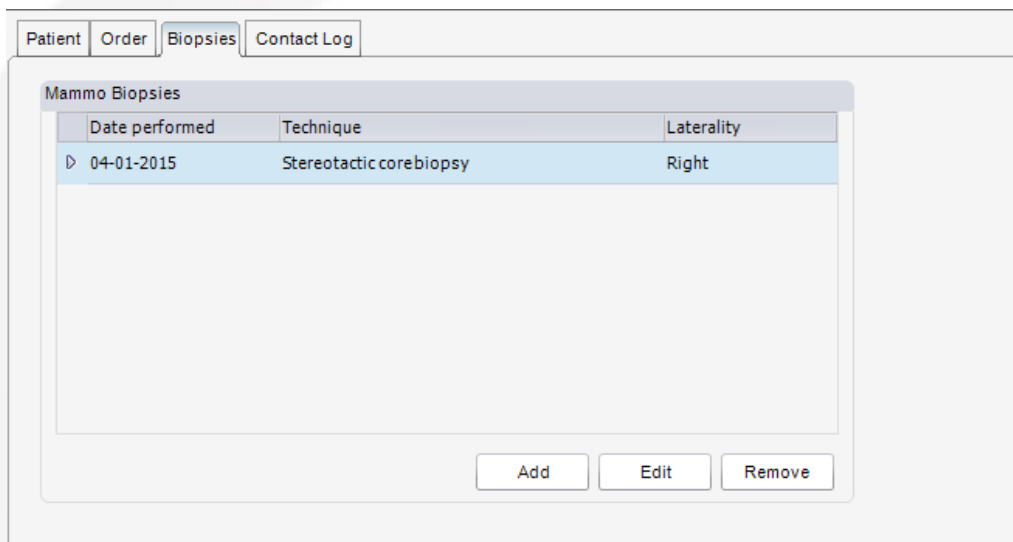

Date performed	Technique	Laterality
04-01-2015	Stereotactic corebiopsy	Right

Then click Save at the bottom of the Follow-Up tab.

- ☞ As long as you do *not* have “Additional Surgery Required” checked off, the exam will drop off of the BI-RADS 4/5 WL.
- ☞ The results will now be included in any MQSA statistics and can be accessed by opening the Patient Folder and opening the exam which recommended the biopsy in View/Edit mode. The results can be found on the Mammography tab.
  - This tab can also be accessed to add pathology results for a breast imaging exam that was not a BI-RADS 4 or 5, if you become aware of a malignancy.
  - When accessed via the Patient Folder, there is an additional option for entering pathology. The Cancer Finding box next to the Mammo Biopsies grid can be used when you do not have enough information to complete the normal biopsy entry, as is sometimes the case with incidental cancers. You do not need to use both.

### Mammo Biopsies

Date performed	Technique	Laterality
----------------	-----------	------------

### Cancer Finding

Positive cancer finding

Date of finding

