


1. Open a new template.

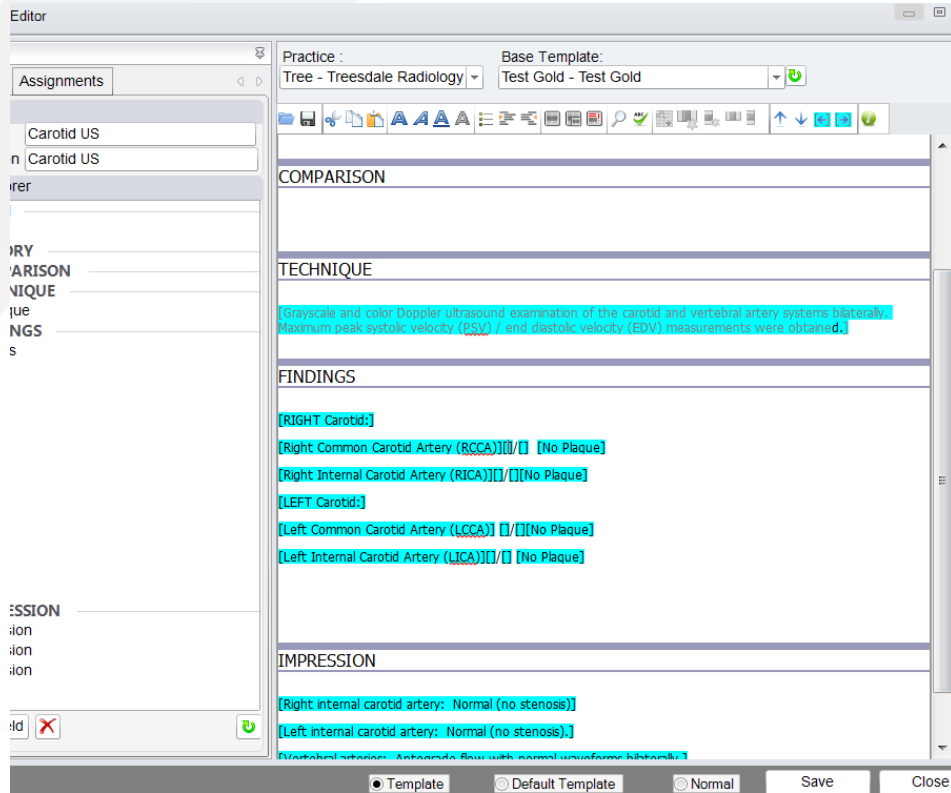

2. Add a field, name it, and click the Advanced button.

The screenshot shows the 'Field Editor' window with the following configuration:

- Name: data6
- Cue: Enter voice navigation cue
- Default: Enter field's default value
- Type: Data field, Text, Pick list, Conditional Val...
- Required field:

Custom field mappings section:

- Currently using: Search fields
- Data fields: Chief complaint, Compare to studies dated..., Contrast administered..., Exam title, Laterality description, Magnet strength, Modality report description, Now (12 hour), Now (24 hour), Patient age in years, Patient gender, Procedure description, Referring name (long), Referring name (short)
- Auto-macros: ExamTitle, Technique
- Advanced field: Dataset, Table, Column, Advanced...

3. Select the Digital Form tab.

The screenshot shows the 'Field Editor' window with the 'Digital Forms' tab selected. The 'Digital Forms' list includes:

- 01 Dave Test Form
- Biopsy Report
- Breast Cancer Risk Assessment
- DC Carotid Test
- Doppler
- eGFR Calc
- FNL Under Mammo/MR
- FNL Under US
- LDCT Scheduling Form
- LDCT Tech
- Mammo Questionnaire
- MammographyQ
- MR/CT Tech Questionnaire
- MRI Screening Questionnaire
- RB - Core Biopsy
- RB - FNA
- RB - FNL left
- RB - FNL right
- RB - Galactography

The 'Questions' section is currently empty. The window has 'OK' and 'Close' buttons at the bottom right.

4. Select the form and field that you want to map and then click OK.

