

Server Update

For eRAD RIS
Version 1.0
Build 40
Update 40.014

Table of Contents

1. New Features	3
PACS Integration Update	3
Reporting – Pending Dictation WL to Omit Studies not Assigned to a Chair	4
Hide Report Contents on Suspended Studies.....	5
2. New Configuration settings and Changes.....	6
3. Defects Resolved	7
4. Package Contents	8
5. Intended Audience	9
6. Who is affected	9
7. eRAD RIS Release Version Numbers	10
8. Installing	12

1. New Features

PACS Integration Update

At times it is possible that the eRAD RIS application will see a delay in the Reporting window when PACS is waiting for the “close session” reply. To mitigate this delay we have introduced a new ris.exe.config application setting “**ShouldCloseViewSessionOnThread**”. The default value is false, meaning the close session happens on the main application thread and we pause for reply, so the reporting window will see a delay loading next. When set to true it closes on a background task and the reporting window will not see the delay. The next study will load in dictate window, although the PACS may still be loading images.

A new ris.exe.config entry has been created with a default value of false.

```
<add key="ShouldCloseViewSessionOnThread" value="false"/>
```

NOTE: With a value of false, the behaviour will be unchanged compared to previous versions.

In the ris.exe.config file, this new key is to be added to the <appSettings> section.

Before:

```
<appSettings>
  <add key="defaultEndPointName" value="NetTcpBinding_IStudyService" />
  <add key="ClientSettingsProvider.ServiceUri" value="" />
  <!--<add key="RowLevelCacheDataSetNamesCSV" value="*" /> -->
</appSettings>
```

After:

```
<appSettings>
  <add key="defaultEndPointName" value="NetTcpBinding_IStudyService" />
  <add key="ClientSettingsProvider.ServiceUri" value="" />
  <add key="LookupTableTimeToLiveInMinutes" value="5"/>
  <add key="ShouldCloseViewSessionOnThread" value="false"/>
  <!--<add key="RowLevelCacheDataSetNamesCSV" value="*" />-->
</appSettings>
```

RADNET, Inc. – eRAD RIS Server Update

Reporting – Pending Dictation WL to Omit Studies not Assigned to a Chair

The Pending Dictation WL content has been modified. Now if a study does not have chair assigned, and the study is not assigned to me (Radiologist user) I will not see the study on my work list. If the study is assigned to me, but is not assigned a chair, it will appear on my Pending Dictation WL. If the study is assigned to the chair I am currently accessing, but is not assigned to a user, it will appear on my work list. Finally, if the study is assigned to the chair I am currently accessing, but is assigned to a different user, it will not appear on my work list. The tables below describe this in better detail. The values in red are what has changed for this build.

If I am sitting in Chair1, then:

Chair	Assigned To	Include in My Pending Dictation WL		
		<= B40.013	>= B40.014	
			Chair WF Enabled	Chair WF Disabled
<null>	<null>	Yes	No	Yes
<null>	Me	Yes	Yes	Yes
<null>	NotMe	No	No	No
Chair1	Me	Yes	Yes	Yes
Chair1	<null>	Yes	Yes	Yes
Chair1	NotMe	No	No	No
Chair2	<null>	No	No	Yes
Chair2	Me	No	Yes	Yes
Chair2	NotMe	No	No	No

If my chair is not set (viewing all chairs), then:

Chair	Assigned To	Include in My Pending Dictation WL		
		<= B40.013	>= B40.014	
			Chair WF Enabled	Chair WF Disabled
<null>	<null>	Yes	No	Yes
<null>	Me	Yes	Yes	Yes
<null>	NotMe	No	No	No
Chair1	Me	Yes	Yes	Yes
Chair1	<null>	Yes	Yes	Yes
Chair1	NotMe	No	No	No
Chair2	<null>	Yes	Yes	Yes
Chair2	Me	Yes	Yes	Yes
Chair2	NotMe	No	No	No

RADNET, Inc. – eRAD RIS Server Update

Hide Report Contents on Suspended Studies

The report template (when modified, this is a customization that is not packaged in the build) will not show the contents of the report for a study that is in suspended status. Looking at the screenshot below, you can see that the report contents are not displayed, but rather a message informing the viewer that the study is suspended. This will also be the case for suspended studies that may be inadvertently sent to print or fax.

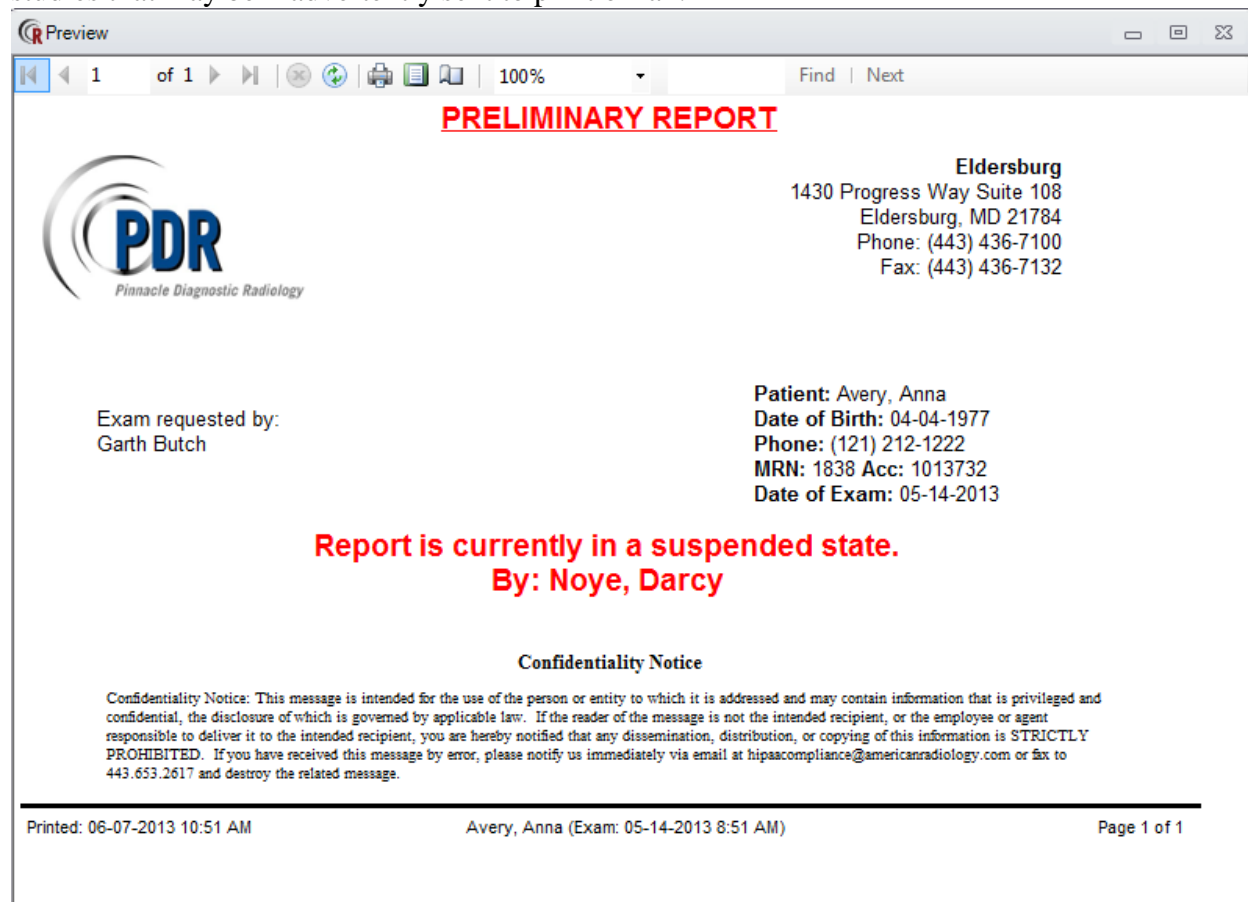


Figure 1.1 – Preview of report in suspended status

Also a new Access String “**Clinical.PreviewSuspendedReport**” has been added to control who has access to preview a report that has been suspended by another user. When a user that does not have this permission tries to preview the report, the option will not be available in the context menu. The data nugget for report will be disabled when the study is in suspended status. Lastly, like the Preview Report context menu option, Send To will not be available, thus removing the option for these users the ability to send the report to Document Distribution.

2. New Configuration settings and Changes

Setting Placeholder	Setting	Default	Purpose / Controls Access to
AccessString lookup UserGroupPermission lookup	Clinical.PreviewSuspended Report	None	Controls access to previewing reports that were suspended by another user
Ris.exe.config	<p>Before:</p> <pre><appSettings> <add key="defaultEndPointName" value="NetTcpBinding_IStudyService" /> <add key="ClientSettingsProvider.ServiceUri" value="" /> <!--<add key="RowLevelCacheDataSetNamesCSV" value="*" /> --> </appSettings></pre> <p>After:</p> <pre><appSettings> <add key="defaultEndPointName" value="NetTcpBinding_IStudyService" /> <add key="ClientSettingsProvider.ServiceUri" value="" /> <add key="LookupTableTimeToLiveInMinutes" value="5"/> <add key="ShouldCloseViewSessionOnThread" value="false"/> <!--<add key="RowLevelCacheDataSetNamesCSV" value="*" /> --> </appSettings> </appSettings></pre>	False	

3. Defects Resolved

Redmine #	Subject	Priority
4019	Patient search error on patient who had multiple merges	Urgent
4091	(Stale cache issue) Can't open scan doc nugget if you have no preference set...	Urgent
3989	Reporting - Error on submit and next	Urgent
3919	RadNet - Prevent users from distributing "SUSPENDED CASES"	Urgent
4040	RadNet - Chair - Don't show null chairs in pending dictation worklist	High
3513	DB showing expensive query on retrieve image of doctors signature for diagnostic report	High
3518	Procedure Plans may not appear correctly on IVT work list	High
4080	logging of failures to obtain a lock	Normal
4075	PACS integration can pause main application thread while waiting for a close session reply	Normal
4069	Memory Leak specific to Reporting and Break Lock	Normal
4064	YellowSticker XML notes sometimes fail and don't properly handle the error	Normal
4059	Reporting - Writing Yellow Sticker not writing exception	Normal
4056	Reporting - Activate/Deactivate errors.	Normal
4050	Reporting - Invoke begininvoke cannot be called...	Normal
4032	Scanning - Canon Flatbed scanners are locking up RIS	Normal
4024	Reporting - Reduce # of WS Calls	Normal
3979	Performance - recently accessed patients should only be saved on app exit	Normal

4. Package Contents

Files included for this server update:

- 1) GUI ZIP Updates folder contains
 - a) rRIS_1.3.40.18063.zip
- 2) RIS Service Updates folder contains
 - a) CDSJobManager.dll
 - b) CDSJobManager.pdb
 - c) CommonLib.dll
 - d) CommonLib.pdb
 - e) DDEngine.dll
 - f) DDEngine.pdb
 - g) rRISServices.exe
 - h) rRISServices.pdb
 - i) rRISServicesx64.exe
 - j) ScheduledTaskManager.dll
 - k) ScheduledTaskManager.pdb
 - l) SignatureApprovalService.exe
 - m) SignatureApprovalService.pdb
 - n) WCFService.dll
 - o) WCFService.pdb
 - p) WebApi.dll
 - q) WebApi.pdb
 - r) Wedge.exe
 - s) Wedge.pdb
- 3) SQL folder contains:
 - a) B40_014.sql

5. Intended Audience

This document is created by the RIS Development team for the RadNet RIS Implementation team.

6. Who is affected

Build 40 installs. This server update must be applied directly to the build 40 core install. Server updates are not cumulative. Any previous server updates must be applied before this one.

7. eRAD RIS Release Version Numbers

Build	Patch	UI Version	Core Version	WS Version	DB Version	Notes
37	-	1.2.37.12843	1.2.37.12782	1.2.37.12843	1.2.37.12843	Full version release
37	1	1.2.37.13040	1.2.37.13040	1.2.37.13039	1.2.37.13003	Core install, GUI.zip, web service and database updates
37	2	1.2.37.13069	1.2.37.13040	1.2.37.13069	1.2.37.13003	GUI.zip and Web service updates
37	3	1.2.37.13069	1.2.37.13040	1.2.37.13069	1.2.37.13003	WebAPI updated, no version changes
37	4	1.2.37.13171	1.2.37.13040	1.2.37.13165	1.2.37.13171	GUI.zip, web service, and database updates
37	5	1.2.37.13347	1.2.37.13040	1.2.37.12979	1.2.37.12979	Bug was discovered with version numbers, this is why it appears we went backwards in the version number but we did not.
37	6	1.2.37.13665	1.2.37.13040	1.2.37.13580	1.2.37.12979	GUI.zip and Web service updates
37	7	1.2.37.13807	1.2.37.13040	1.2.37.13799	1.2.37.13807	GUI.zip and Web service updates. Also MIRTH mapping changes
37	8	1.2.37.13807	1.2.37.13040	1.2.37.13827	1.2.37.13807	Web service update
37	9	1.2.37.13807	1.2.37.13040	1.2.37.13839	1.2.37.13807	Web service update
37	10	1.2.37.13895	1.2.37.13040	1.2.37.13895	1.2.37.13895	GUI.zip, Web service and database updates
37	11	1.2.37.14161	1.2.37.13040	1.2.37.14161	1.2.37.14161	GUI.zip, Web service, database updates and SQL Management reporting
37	12	1.2.37.14267	1.2.37.13040	1.2.37.14161	1.2.37.14267	GUI.zip, Web service and database updates
37	13	1.2.37.14326	1.2.37.13040	1.2.37.14326	1.2.37.14326	GUI.zip, Web service and database updates
38	-	1.2.38.14486	1.2.38.14486	1.2.38.14486	1.2.38.14486	Full version release
38	1	1.2.38.14748	1.2.38.14748	1.2.38.14748	1.2.38.14748	Core install required
38	2	1.2.38.15182	1.2.38.14748	1.2.38.14748	1.2.38.14748	GUI.zip update only
38	3	1.2.38.15251	1.2.38.14748	1.2.38.14748	1.2.38.14748	GUI.zip update only
38	4	1.2.38.15327	1.2.38.14748	1.2.38.15290	1.2.38.14748	GUI.zip and Web service updates
38	5	1.2.38.15513	1.2.38.14748	1.2.38.15513	1.2.38.15451	GUI.zip, Web service and database updates
38	6	1.2.38.15618	1.2.38.14748	1.2.38.15513	1.2.38.15618	GUI.zip and database updates
38	7	1.2.38.15768	1.2.38.14748	1.2.38.15768	1.2.38.15618	GUI.zip and Web Service updates
39	-	1.3.39.15389	1.3.39.15389	1.3.39.15377	1.3.39.15355	Full version release
39	1	1.3.39.16026	1.3.39.15943	1.3.39.15943	1.3.39.15957	Requires Core Install. GUI.zip, Web Service, and Database updates
40	-	1.3.40.16628	1.3.40.16628	1.3.40.16628	1.3.40.16628	Full version release
40	1	1.3.40.16849	1.3.40.16659	1.3.40.16744	1.3.40.16727	GUI.zip, Web service, database updates and SQL Management reporting
40	2	1.3.40.16868	1.3.40.16659	1.3.40.16868	1.3.40.16727	GUI.zip and Web Service

RADNET, Inc. – eRAD RIS Server Update

40	3	1.3.40.16868	1.3.40.16659	1.3.40.16868	1.3.40.16727	DB updated only
40	4	1.3.40.16967	1.3.40.16659	1.3.40.16932	1.3.40.16961	GUI.zip, Web service, database updates
40	5	1.3.40.17087	1.3.40.16659	1.3.40.17069	1.3.40.17087	GUI.zip, Web service, database updates
40	6	1.3.40.17146	1.3.40.16659	1.3.40.17146	1.3.40.17146	GUI.zip, Web service, database updates
40	7	1.3.40.17268	1.3.40.16659	1.3.40.17268	1.3.40.17259	GUI.zip, Web service, database updates
40	8	1.3.40.17364	1.3.40.16659	1.3.40.17362	1.3.40.17343	GUI.zip, Web service, database updates
40	9	1.3.40.17370	1.3.40.16659	1.3.40.17370	1.3.40.17370	GUI.zip, Web service, database updates
40	10	1.3.40.17435	1.3.40.16659	1.3.40.17435	1.3.40.17370	64 BIT web service. GUI.zip, Webservice.
40	11	1.3.40.17435	1.3.40.16659	1.3.40.17525	1.3.40.17370	Web service update
40	12	1.3.40.17685	1.3.40.16659	1.3.40.17685	1.3.40.17648	GUI.zip, Web service, database updates
40	13	1.3.40.17776	1.3.40.16659	1.3.40.17776	1.3.40.17648	GUI.zip, Web service
40	14	1.3.40.18063(3GB)	1.3.40.16659	1.3.40.18063	1.3.40.18063	GUI.zip, Web service, database updates

8. Installing

- 1) Copy and replace the current rRIS...zip file with the rRIS_1.3.40.18063.zip file provided with this release.

Be sure to make a backup of the ris.exe.config file contained within the current .zip folder.

NOTE: there is a change to the ris.exe.config file.

See Sections

1 New Features - Configurable dataset retrieval interval and

2 New Configuration settings and Changes

- 2) Stop each service listed below. Copy the files from the RIS Service Updates folder to the directory for each rRISService. Be sure to make backup of all files, but especially the rRISServices.exe.config and rRISServices.vshost.exe.config files.

rRISService directories include:

- a) *rRISService*
- b) *CDS rRISService*
- c) *Wedge rRISService*
- d) *DDE rRISService*

Restart each service.

- 3) Run the B40_014.sql file against the eRAD RIS database instance.